


SYNCA

## User Manual

Quzy

Neck & Shoulder Massager

Model: SM260

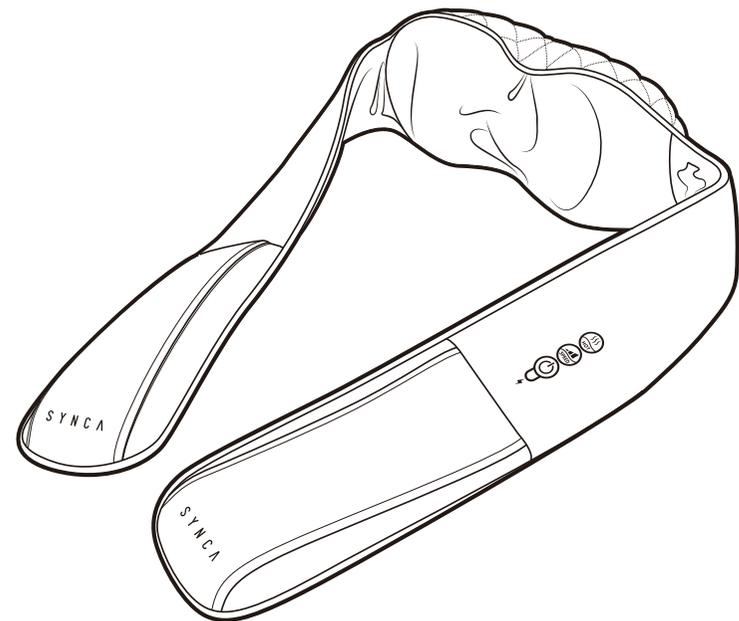

This manual should be kept on hand for easy reference.  
It is your guide to a safe and efficient operation.

## IMPORTANT SAFETY INSTRUCTIONS

When using the product, all instructions and basic precautions should always be followed, including the following:

### DANGER

– To reduce the risk of electric shock:  
Always unplug this product from the electrical outlet immediately after use and before cleaning.

### WARNING

- To reduce the risk of burns, fire, electric shock, or injury to persons:
  - 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not use and before putting on or taking off parts.
  - 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
  - 3) Close supervision is necessary when this appliance is used by, or near children, invalids, or disabled persons.
  - 4) Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
  - 5) Never operate this product if it is dropped, has a damaged cord or plug, a malfunction, or it is dropped into water. Return the product to service center for examination and repair.
  - 6) Do not carry this product by its supply cord or use cord as a handle.
  - 7) Keep the cord away from heated surfaces.

## WARNING

- 8) Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 9) Never drop or insert any object into any opening.
- 10) Do not use outdoors.
- 11) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 12) To disconnect, switch off all controls and remove plug from outlet.
- 13) This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- 14) Children should be supervised to ensure that they do not misuse the product.
- 15) The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- 16) The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- 17) The appliance is only to be used with the power supply unit provided with the appliance.

### SAVE THESE INSTRUCTIONS

This product is intended for household use only

## HOW TO OPERATE



1. ON or OFF and Kneading direction control  
Press to start working, automatically changing kneading direction every minutes, press again to set only one direction, then another direction by the third time. the unit will be shut-off automatically in 15 minutes, Long press the button for shut-off.

Charging indicator:

when plug the adaptor, the light turn red means products under charging, once indicator light turn green means fully charged. If start massage, products will stop charging then the indicator light goes out.




2. Speed control  
3 levels of speed, the default is high speed, press once lower speed, press twice middle speed, press again high speed.




3. Heater  
You can choose to stop heating function when it is working, press one time, stop heating function, press again, heating function work again.

4. Press the power button to stop the massage at any time.

## MASSAGE POSITION

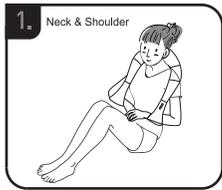


## NAME AND FUNCTION OF THE MACHINE PARTS

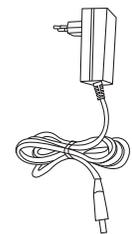

Power adaptor

Rechargeable type

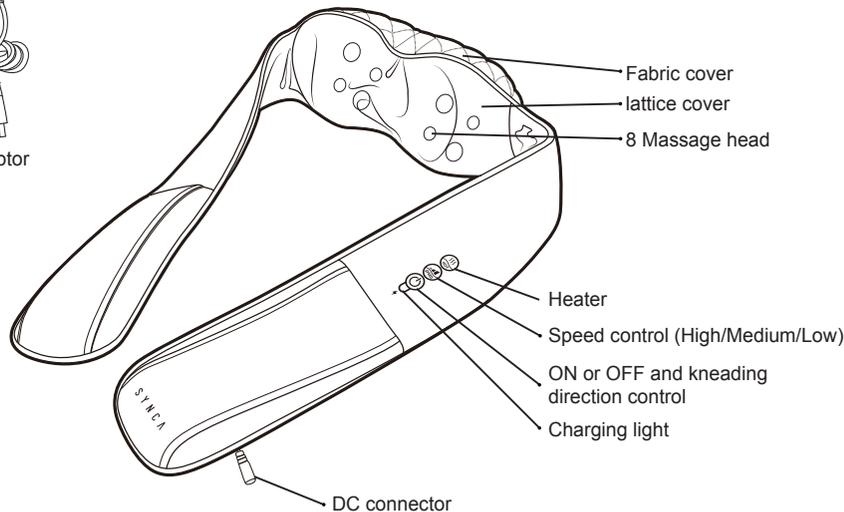

## TO CLEAN AND MAINTAIN THE PRODUCT

### 1. Cleaning

Before you clean, please make sure that the switch is in the off location. And the power plug has been pulled out from the socket. Please use dry cloth to clean the product.

### ⚠ Attention

- Please make sure that do not use benzene or thinner and other aggressive detergents to clean the product. Or the product may fade.

### 2. Maintenance

Make sure that the switch is off, and the power plug has been pulled out from the socket. Keep the product clean and far from the wet environment.

## TECHNOLOGY PARAMETER

- Name: Neck & Shoulder Massager  
Type: SM260  
Adapter: Input: 100-240V~50/60Hz  
Output: DC 13V 1.8A  
Battery Type: Li-ion 2200mAh 11.1V  
Battery Life: Approx. 2.5 hours (without heat)  
Battery Charge Time: Approx. 2.5 hours